



## I.V. POLE STANDS AND TOOL-FREE ACCESSORIES



Medline has an I.V. Pole Stand to meet all your needs and budgets — from standard duty to fully-featured heavy-duty models.

# PREMIER HEAVY-DUTY I.V. POLE

WITH QUICK RELEASE CASTERS

**Ideal for triple-channel pump or any application that would require an extra stable base**

- Patent-pending 3" (7.6 cm) premium Quick Release Casters help keep pole in service longer
- Heavy-duty 24" (61 cm) diameter low center of gravity hammer-tone cast iron base with 6 wheels help resist tipping, weighs 19.7 lbs. (9 kg)
- Rubber bumpers at end of legs help to protect walls and other objects from scuffing
- D-tube insert helps ensure stability of pole without the fear of base spinning
- Stainless steel telescoping pole adjusts from 48" (122 cm) compressed to 82" (208 cm) fully extended
- Four durable ram hooks

## TOOL-FREE ACCESSORIES FOR MDS80600 MODEL

### URINE BAG HOOK

Item#: MDS80600UBH\*  
Pkg: 10/cs



### STAND ASSIST WITH TRAY

Item#: MDS80600SAT\*  
Pkg: 6/cs



### GRAB BAR

Item#: MDS80600GB\*  
Pkg: 6/cs



### BASKET

Item#: MDS80600B\*  
Pkg: 2/cs



### OXYGEN TANK HOLDER

Item#: MDS80600OTH\*  
Pkg: 2/cs



### CASTERS

Item#: MDS80600C\*  
Pkg: 2/bx



ADD VERSATILITY

\* Patent Pending

ONE OF THE  
MOST STABLE BASES  
AVAILABLE!



*pump is not included*

MDS80600SAT  
(optional)



MDS80600B  
(optional)

MDS80600UBH  
(optional)

ITEM# MDS80600  
PKG. 2/CS



TO ORDER: 1-800-MEDLINE


## HEAVY-DUTY ALUMINUM

### I.V. POLE STAND

**Ideal for pumps, multiple bags, or any application which requires an extra stable base**

- Constructed of corrosion-resistant aluminum with heavy-duty 18 lb. (8.2 kg) cast iron five-legged base with low center of gravity to help ensure stability
- Base diameter 21.5" (54.6 cm)
- Accommodates two pumps up to 14 lbs. (6.3 kg) each and two 3-liter bags of fluid for versatility
- Friction grip lock makes adjusting easy
- Height adjustable from 46.75" to 81" (119 cm to 205 cm)
- Four durable ram hooks hold items securely

ITEM# MDS80500  
PKG. 2/CS

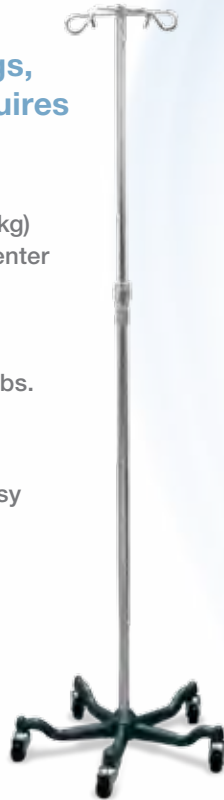


## DELUXE

### I.V. POLE STAND

**Standard duty model featuring an extra-high height range**

- Constructed of corrosion-resistant aluminum with five-legged, 23.5" diameter (59.7 cm) cast aluminum base weighing 10 lbs. (4.5 kg)
- 2" (5 cm) heavy-duty casters provide smooth, quiet movement
- Accommodates two pumps up to 14 lbs. (6.3 kg) and two 3-liter bags of fluid for versatility
- Microphone-style knob makes adjustments easy
- Height adjustable from 56.75" to 96" (144 cm to 244 cm)
- Four ram hooks hold items securely

ITEM# MDS80494  
PKG. 2/CS

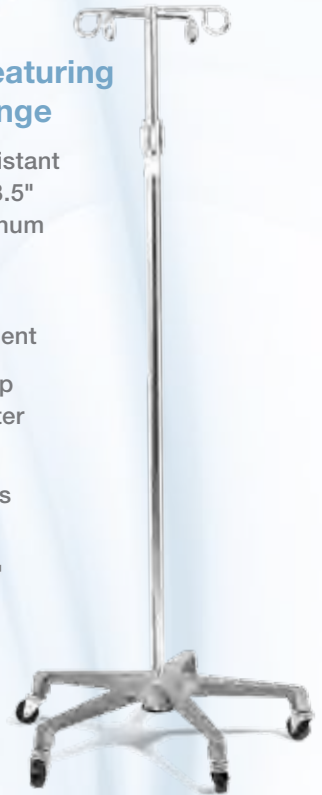

## STANDARD

### I.V. POLE STAND

**Value-positioned to meet basic needs**

- Large 23" (58 cm) diameter, four-legged base provides excellent stability
- Height adjustable from 47" to 85" (119.4 cm to 216 cm)
- Base and pole made of stainless steel for durability
- Two ram hooks hold solution containers securely

ITEM# MDS80441  
PKG. 1 EACH



## HOMECARE

### I.V. POLE STAND

**Economical for light-duty home use**

- Five-legged plastic base offers greater stability
- Base diameter of 23.5" (59.7 cm)
- Tool-free push-button design makes adjusting easy
- Height adjustable from 37.5" to 83.5" (95.3 cm to 212 cm)
- Two ram hooks hold items securely

ITEM# MDS129CS  
PKG. 3/CS



## OTHER POPULAR PRODUCTS

# STAINLESS STEEL

ITEM# MPH08MS01C  
PKG. 1 EACH



### Mayo stand

- Tray measures 16" x 21" (41 cm x 53 cm)
- Foot operated for ease and convenience
- Height adjusts from 39.5" to 62" (100 cm to 157 cm)

ITEM# MPH08190  
PKG. 1 EACH



### Kickbucket

- 13 qt (12.3 liter)
- With wall guard rubber bumper and 4 roller blade wheels for smooth, wobble-free movement

ITEM# MPH08SS01C  
PKG. 1 EACH



### Stacking stool

- 12"L x 18"D x 6"H (30 cm x 46 cm x 15 cm)
- Stability bars for extra safe footing
- Connection points for secure assembly of multiple units

ITEM# MPH0807SSHB  
PKG. 1 EACH



### Single basin stand with H-brace

- 8.5 qt (8 liter) basin
- 15"L x 15.25"W x 33"H (38 cm x 39 cm x 84 cm)

## GET ALL THE DETAILS

Contact your Medline sales representative to find out how we can help you with your equipment needs.



**Medline Industries, Inc.**  
One Medline Place  
Mundelein, IL 60060

Medline United States  
**1-800-MEDLINE (633-5463)**

[www.medline.com](http://www.medline.com) | [info@medline.com](mailto:info@medline.com)

Medline Canada  
**1-800-396-6996**

[www.medline.ca](http://www.medline.ca) | [canada@medline.com](mailto:canada@medline.com)

Medline México  
**01-800-831-0898**

[www.medlinemexico.com](http://www.medlinemexico.com) | [mexico@medline.com](mailto:mexico@medline.com)