

INVACARE®

# microAIR® MA50

SERIES THERAPEUTIC SUPPORT SURFACES

---



The **Invacare® microAIR® MA50 Therapeutic Support Surface** overlay is a full-featured alternating pressure overlay with four alternating pressure cycles to provide relief and comfort.

The **Invacare® microAIR® MA51 Therapeutic Support Surface** model is an alternating pressure mattress replacement system which helps reduce interface pressure to aid in the prevention and treatment of Stage I to IV pressure ulcers.

In addition to all of the features of the microAIR MA51 mattress, the **Invacare® microAIR® MA55 Therapeutic Support Surface** model improves air circulation and cools the body temperature with a unique system of low air loss released into the mattress.

## Features

**State of the art compressor** – With solenoid valve system allows for easy set up and effective therapy cycles  
Four therapy times – Can be adjusted for patient comfort to provide pressure relief

**Constant pressure monitoring** – Constantly monitors the air pressure in the mattress to maintain the appropriate pressure to prevent capillary occlusion (generally below 32mm/Hg)

**Manual Fowler** – When activated manually, 25% more air is automatically pushed into the mattress to prevent patient from bottoming out when bed is articulated

**Audio and visual alarms** – Warns the patient and caregiver during power outage or low pressure



**Lock out button** – prevents against tampering or accidental adjustments

**Quick connect coupler** – allows for quick and easy set up

**Power outage safety** – Mattress remains inflated during power outage as long as air hose is connected to the compressor

**Pressure adjustment** – Pressure can be adjusted by 25%, 50% or 75% to provide comfort for frail patients

**HCPSC Code: E0372(MA50), E0277(MA51/MA55)**

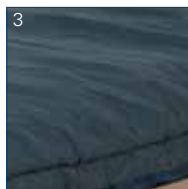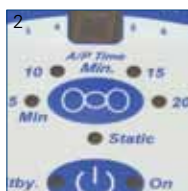

1. Easy to replace cells in case of damage
2. Choice of four therapy cycles
3. Dual layer quilted cover, denier nylon bottom



**Invacare Corporation**  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2012

Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
Form No. 10-125  
Rev. 03/12  
111232



**Yes, you can.®**

Model #	Dimensions	Weight Capacity	Compressor	Alternating Time Cycle	Alarms	Limited Warranty
MA50	80" L 36" W 6" H	350 lb.	10 LPM	Adjustable	Audio and Visual	1 Year (Cover/Mattress) 2 Years (Power Unit)
MA51	80" L 36" W 8" H	350 lb.	10 LPM	Adjustable	Audio and Visual	1 Year (Cover/Mattress) 2 Years (Power Unit)
MA55	80" L 36" W 8" H	350 lb.	10 LPM	Adjustable	Audio and Visual	1 Year (Cover/Mattress) 2 Years (Power Unit)

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.