

INVACARE®

microAIR® MA90

SERIES THERAPEUTIC SUPPORT SURFACE



The **Invacare® microAIR® MA90Z Therapeutic Support Surface** model is an advanced mattress enhanced with patented features to create the ultimate patient comfort and relief. Using a compressor system, the MA90Z mattress offers alternating lateral rotation and on demand low air loss to help provide pressure relief and moisture management for up to stage IV pressure ulcers.

The **Invacare® microAIR® MA95Z Therapeutic Support Surface** model allows customization for degrees of turn, time and pressure. Static and pulse therapies are available in addition to the alternating lateral rotation and continuous true low air loss. Added to the MA90Z and MA95Z are enhanced protective measures such as side bolsters and recentering of the patient in between turns. The MA95Z aids in the prevention and treatment of stage I to IV pressure ulcers.

Features

MA90Z Compressor – State of the art compressor allows for easy set up and effective therapy cycles

MA90Z Blower – Powerful blower provides greater air output, quick set up and more efficient therapy cycles

Constant pressure monitoring – Solenoid valves open and close to regulate the air inside the mattress and maintain the pressure at or below 32mm/Hg to avoid capillary occlusion

Wireless auto fowler – When the transmitter in the head section is raised, 25% more air is automatically added into the mattress, helping prevent the patient from bottoming out


Power outage safety – Compressor model (MA95Z) remains inflated through power outage

Turn times – MA90Z includes times of 10, 20, 30 and 60 minutes with option to turn right or left only. MA95Z includes customizable turn times

Rotation angles – MA90Z includes ¼, ½ and ¾ turn options. MA95Z features customizable turn angles

Z-Cell – Patented design creates a more comfortable and subtle turn

HCPCS Code: E0277

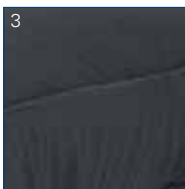

1. 24 air cells and built-in raised side bolsters
2. Audio and visual alarms with lock out
3. Urethane four way stretch quilted cover, denier nylon bottom



Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2012

Invacare Corporation.
All rights reserved.
Trademarks are identified
by the symbols ™ and
®. All trademarks are
owned by or licensed
to Invacare Corporation
unless otherwise noted.
Specifications are subject
to change without
notification.
Form No. 10-132
Rev. 03/12
111235



Model #	Dimensions	Weight Capacity	Compressor	Blower	Alarms	Regulation	Limited Warranty
MA90Z	80" L 36" W 10" H	350 lb.	50 LPM	N/A	Audio and Visual	16 CFR 1632 16 CFR 1633	1 Year(Cover/ Mattress) 2 Years (Power Unit)
MA90ZB42/ MA90ZB48	80" L 42"/48" W 10" H	1,000 lb.	50 LPM	N/A	Audio and Visual	16 CFR 1632 16 CFR 1633	1 Year(Cover/ Mattress) 2 Years (Power Unit)
MA95Z	80" L 36" W 10" H	350 lb.	N/A	1275 LPM	Audio and Visual	16 CFR 1632 16 CFR 1633	1 Year(Cover/ Mattress) 2 Years (Power Unit)
MA95ZB42/ MA95ZB48	80" L 42"/48" W 10" H	1,000 lb.	50 LPM	N/A	Audio and Visual	16 CFR 1632 16 CFR 1633	1 Year(Cover/ Mattress) 2 Years (Power Unit)

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.