

BLICKMAN



Decontamination Stations

Includes:

- Decontamination Stations
- Faucet Options
- Accessories

Optimal infection control

OUR MISSION

Customer Focus

Exceed customer expectations through a commitment to continuous improvement of our products and services.

Quality Products

Offer the marketplace a broad range of high quality, cost effective, and durable products of superior craftsmanship.

Organizational Excellence

Strive for excellence in all we do by fostering a participative, team based environment, which nurtures the active involvement of all members of the organization.

Marketplace Recognition

Earn the respect, admiration, and loyalty of our people, customers, and suppliers by treating them with dignity and respect, embracing diversity, and developing and maintaining long term partnerships.

Collaborative Workplace

Encourage all employees to work safely and efficiently, contribute innovative ideas, and personally benefit from the ongoing success of the company.


The logo for Blickman, featuring the word "BLICKMAN" in a bold, white, sans-serif font. The text is contained within a dark blue rectangular box with a thin white border. The box has a slight 3D effect with a lighter blue shadow on the right side.

THE BLICKMAN DIFFERENCE

Decontamination Stations

- Our decontamination stations provide optimal infection control in today's high risk healthcare environments. Our durable stainless steel decontamination stations meet the stringent AAMI (Association for the Advancement of Medical Instrumentation) guidelines for maximum sterility and utility.
- Our decontamination stations have an ergonomic working height, symmetrical design for either left or right work flow, 18" counter space on both sides, and an 8" backsplash.
- Three sink bowls are included, which feature a unique, "no drip" sink lip, pitched bottoms for proper drainage, and permanently etched gallon and liter calibration inside each bowl for precise measurement.
- Faucets include: a pre-rinsed combination faucet, a wall mounted faucet, and a wall mounted air, gas and VAC valve with serrated barbed end hose.

Key Product Features

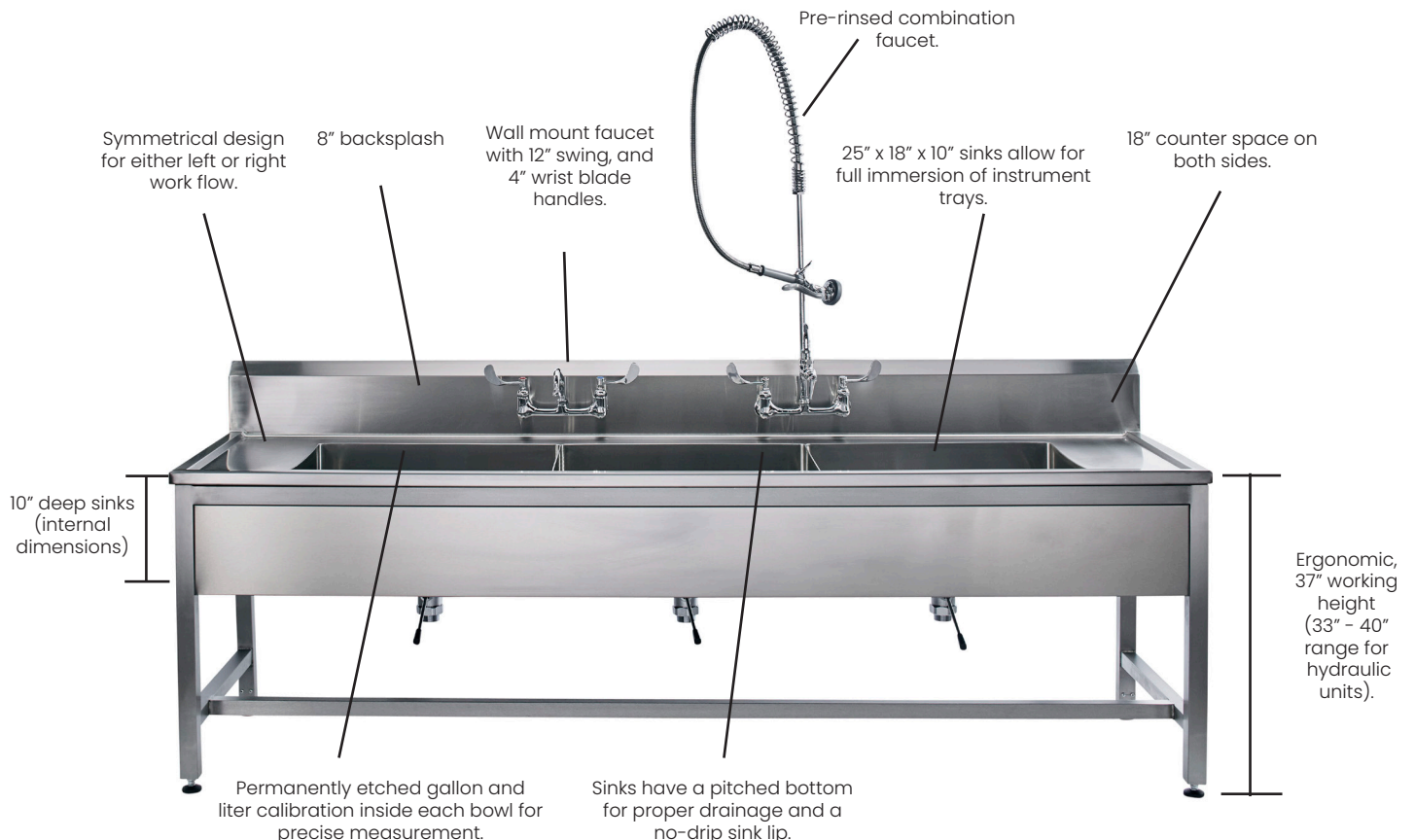

Decontamination Stations

Description: We offer three decontamination stations. One sits on two base cabinets, and the other two sit on stainless steel legs that are either stationary or adjustable. All three types of decontamination stations are set to ergonomic working heights. You have the choice of wall mounted and deck mounted faucets.

Decontamination Station External Dimensions

Width: 116" x Height: 37" x Depth: 26 1/2" or 30"

*33" - 40" adjustable height range for hydraulic unit



Decontamination Stations

Two Base Cabinets

Mounting Types	Wall mounted Faucets 26 1/2" Deep	Deck Mounted Faucets 30" Deep
Model #	DSBC26	DSBC30
Item #	162DSBC001	162DSBC000
Shipping Weight	371 lbs.	391 lbs.

Stationary Legs

Mounting Types	Wall mounted Faucets 26 1/2" Deep	Deck Mounted Faucets 30" Deep
Model #	DSHB26	DSHB30
Item #	162DSHB001	162DSHB000
Shipping Weight	442 lbs.	462 lbs.

Hydraulic Legs

Mounting Types	Wall mounted Faucets 26 1/2" Deep	Deck Mounted Faucets 30" Deep
Model #	DSHY26	DSHY30
Item #	162DSHY017	162DSHY000
Shipping Weight	442 lbs.	462 lbs.

All decontamination stations include three station sinks.

Sink Dimensions: 25" x 18" x 10".



Decontamination Stations

Options

FAUCET OPTIONS

Deck Mounted



Wall Mounted



Name	Pre-rinse Combo Faucet	Swing Faucet	Swivel Gooseneck Faucet	DI/RO Lab Faucet & Lab Hose Fitting Kit
Deck Mount				
Part #	9200601000	9200002214	9243232204	9200586002
Shipping Weight	13 lbs.	5 lbs.	5 lbs.	5 lbs.
Wall Mount				
Part #	9200133011	9200231021	9243232414	9204250KIT
Shipping Weight	13 lbs.	5 lbs.	5 lbs.	2 lbs.

Note: The DI/RO Lab Faucet is not for water, but for gas, air, or vacuum hose purposes.



Decontamination Stations

Accessories



Peg Board
(Image is representative)



Eyewash



Overflow Back Drain




Cleaning Spray Gun



Thermo Valve

Name	Peg Board (24" x 116")	Eyewash	Overflow Back Drain	Cleaning Spray Gun + Tabletop Rosette	Thermo Valve
Part #	2701160240	9201849LHL	9201074150	9200107005 + 9200800017	9203600LF0
Shipping Weight	143 lbs.	8 lbs.	6 lbs.	2 lbs. + 1 lb.	7 lbs.

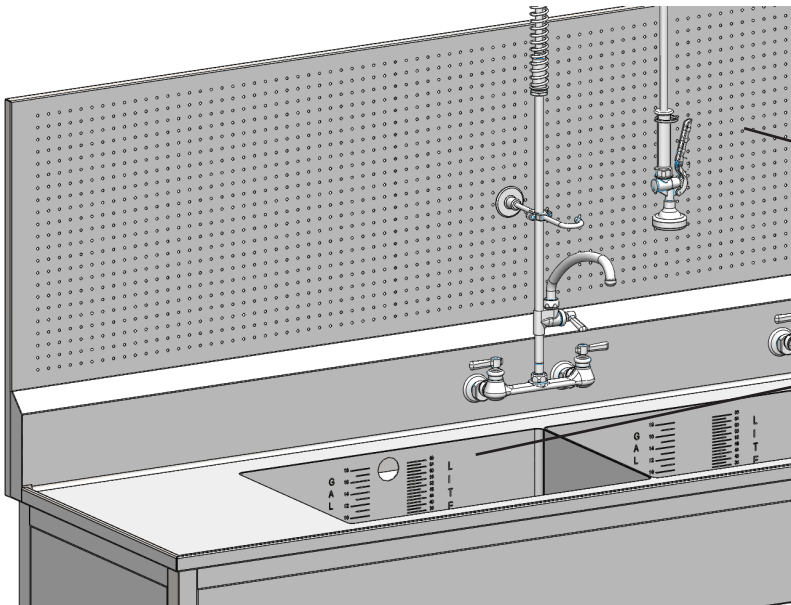

Accurate Diagram of Full Pegboard Back.

Etched Measurements in Sink Bowl

Note: The Cleaning Spray Gun comes with 1 tip. 8 different sized replacement tips can be ordered:
Part# 9200107006.



Decontamination Stations

Construction Specifications

Construction Specifications

Standard Dimensions	External Dimensions: Width: 116" x Depth: 26 1/2 or 30" x Height: 37" (33" - 40" hydraulic unit) Internal Sink Bowl Dimensions: Width: 25" x Height: 18" x Depth: 10"
Materials	304 stainless steel
Gauge	Bowl: 14GA Sink Top and Channel: 16GA DSHY Legs: 16GA; Cross Bracing Square Tube: 18GA DSHB Legs and Cross Bracing Square Tube
Welding	All welding by heliarc method and exposed welds are ground smooth.
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish.
Electrical	Hydraulic Unit: 120V AC.
Casters	N/A
Notes	Ergonomic working height, symmetrical design for either left or right work flow, at least 18" counter space on both sides, and an 8" backsplash. Three sink bowls which feature a unique, "no drip" sink lip, pitched bottoms for proper drainage, and permanently etched gallon and liter calibration inside each bowl for precise measurement.

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.



OUR PROMISE

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality, and we pride ourselves on providing top notch and personalized service to ensure our customers are 100% satisfied. If you are ever not completely satisfied for any reason, we will make it right. And with our equipment backed by our Limited Lifetime Warranty, you can be certain you have made the very best choice for your healthcare facility.