

BLICKMAN



Exam Stools

Providing ultimate comfort and function

Includes:

- Spindle Adjusted Stools
- Pneumatic Stools
- Parts & Accessories

OUR MISSION

Customer Focus

Exceed customer expectations through a commitment to continuous improvement of our products and services.

Quality Products

Offer the marketplace a broad range of high quality, cost effective, and durable products of superior craftsmanship.

Organizational Excellence

Strive for excellence in all we do by fostering a participative, team based environment, which nurtures the active involvement of all members of the organization.

Marketplace Recognition

Earn the respect, admiration, and loyalty of our people, customers, and suppliers by treating them with dignity and respect, embracing diversity, and developing and maintaining long term partnerships.

Collaborative Workplace

Encourage all employees to work safely and efficiently, contribute innovative ideas, and personally benefit from the ongoing success of the company.

The logo for Blickman, featuring the word "BLICKMAN" in a bold, white, sans-serif font. The text is contained within a dark blue rectangular box with a thin white border. The box is set against a light gray background that has a subtle gradient and a slight shadow effect.

THE BLICKMAN DIFFERENCE

- Our durable, high quality chrome and stainless steel exam stools come with 3" thick, dense foam 15" diameter cushion for ultimate comfort.
- Our upholstery is washable Vinyl and is available in an array of 74 fashion colors that meet specific design requirements.
- For applications that require the longevity that stainless steel provides, pick from our revolving stools, or for a more modern look, choose from four or five leg models with spindle adjusted or pneumatic lift mechanisms.
- For back support, choose one of our models with a back rest.
- Our three styles of locking casters allow for easy mobility and safety.




Exam Stools

Spindle

4 LEG SPINDLE EXAM STOOLS

Description: We offer two four leg spindle adjusted exam stools with the option of 2" casters or rubber feet tips. Choose our passaic stool and sit comfortably on a 15" diameter, 3" thick polyfoam padded seat upholstered in vinyl or choose our Clifton stool for an easy-to-clean stainless steel seat with a non-slip vinyl insert. Add a removable pad slipcover to your Clifton stool for a more comfortable option.

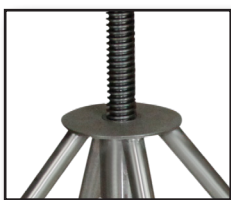

Spindle

Adjusts 12"

Rubber Tips: 19" - 31"

Casters: 21" - 33"



| Name | CLIFTON | | PASSAIC |
|-----------------|------------|------------|------------|
| | Standard | With Pad | Standard |
| | Casters | | |
| Model # | 7745SS-C | 7745SS-PC | 7714-PSSC |
| Item # | 1027745001 | 1027745003 | 1027714001 |
| Shipping Weight | 17 lbs. | 17 lbs. | 17 lbs. |
| | Feet | | |
| Model # | 7745SS | 7745SS-P | 7714-PSS |
| Item # | 1027745000 | 1027745002 | 1027714000 |
| Shipping Weight | 16 lbs. | 16 lbs. | 16 lbs. |

Options: Available in standard black. Consult factory for requests in other colors for large order quantities.



Exam Stools

Spindle

5 LEG SPINDLE EXAM STOOL

Description: Our five leg spindle exam stool rests on a 24" diameter base. Choose between a black composite base or an aluminum base with Wall Saver™ tips on the ends to protect against scuff marks and offer a textured surface for foot grip. Sit comfortably on its 15" diameter, vinyl upholstered seat with 3" thick polyfoam padding. This model is available with or without a back.



Spindle

Adjusts 5 7/8"
Height: 20 7/8" - 26 3/4"



| Name | 5 Leg Spindle Stool | | 5 Leg Spindle Stool with backrest |
|-----------------|----------------------|--|-----------------------------------|
| | Aluminum Base | | |
| Model # | 1113 | | 1113WB |
| Item # | 1021113000 | | 1021113200 |
| Shipping Weight | 22 lbs. | | 26 lbs. |
| | Black Composite Base | | |
| | | | |
| Model # | 1113BLK | | 1113WB-BLK |
| Item # | 1021113100 | | 1021113300 |
| Shipping Weight | 22 lbs. | | 26 lbs. |

Options: Foot ring, 74 different Vinyl colors.



Exam Stools

Pneumatic

FOOT OPERATED EXAM STOOL

Description: Our five leg foot operated exam stool have a push ring under a 24" diameter aluminum base with Wall Saver™ tips on the ends to protect against scuff marks and offer a textured surface for foot grip. Sit comfortably on its 15" diameter, vinyl upholstered seat with 3" thick polyfoam padding. This model is available with or without a back.



Foot Pedal

Adjusts 7"
Height: 18 3/8" - 25 3/8"



| Name | Foot Operated Stool | Foot Operated Stool with Backrest |
|-----------------|---------------------|-----------------------------------|
| Aluminum Base | | |
| Model # | 1205 | 1206 |
| Item # | 1041205000 | 1041206000 |
| Shipping Weight | 22 lbs. | 28 lbs. |

Options: Foot ring, 74 different Vinyl colors.



Exam Stools

Pneumatic

HAND OPERATED EXAM STOOLS

Description: Our five leg hand operated exam stool rests on a 24" diameter base. Choose between a black composite base or an aluminum base with Wall Saver™ tips on the ends to protect against scuff marks and offer a textured surface for foot grip. Sit comfortably on its 15" diameter, vinyl upholstered seat with 3" thick polyfoam padding. This model is available with or without a back.

Standard Stools



Hand Lift

Adjusts 7"
Height: 18 3/8" - 25 3/8"



| Name | Hand Operated Stool | Hand Operated Stool with Backrest |
|-----------------|----------------------|-----------------------------------|
| | Aluminum Base | |
| Model # | 1210 | 1210WB |
| Item # | 1041210000 | 1041210200 |
| Shipping Weight | 22 lbs. | 26 lbs. |
| Name | Black Composite Base | |
| | Hand Operated Stool | Hand Operated Stool with Backrest |
| Model # | 1210BLK | 1210WB-BLK |
| Item # | 1041210100 | 1041210300 |
| Shipping Weight | 22 lbs. | 26 lbs. |

Options: Foot ring, 74 different Vinyl colors.



Exam Stools

Pneumatic

HAND OPERATED EXAM STOOLS

Economy Stools

Description: Our economy hand operated exam stool rests on a 24" diameter black composite base with black plastic twin wheel casters.. Sit in comfort on a 15" diameter, vinyl upholstered seat with 3" thick polyfoam padding. this model is available with and without a back.

Adjusts 7" from 18 ³/₈" to 25 ³/₈" off the ground.



| Name | Economy Stool | |
|-----------------|---------------|------------|
| | Without Back | With Back |
| Model # | 1200 | 1201 |
| Item # | 1041200100 | 1041201300 |
| Shipping Weight | 22 lbs. | 26 lbs. |

Options: Foot ring, 74 different Vinyl colors.

Deluxe Stools

Description: Our deluxe hand operated exam stool comes standard with a foot ring for leg-joint comfort. These stools rest on a 24" diameter black composite or aluminum base, which comes with Wall Saver™ tips to protect against scuff marks and offer a textured surface for foot grip. Both bases travel easily on 2" hard rubber casters. Sit in comfort on a 15" diameter, vinyl upholstered seat with 3" thick polyfoam padding with a backrest.

Adjusts 7" from 21 ³/₈" to 28 ³/₈" off the ground.



| Name | Deluxe Stool | |
|-----------------|---------------|----------------------|
| | Aluminum Base | Black Composite Base |
| Model # | 1212 | 1212BLK |
| Item # | 1041212000 | 1041212100 |
| Shipping Weight | 29 lbs. | 29 lbs. |

Options: 74 different Vinyl colors.

Foot Ring as an Accessory: Compatible on all 5 leg stools

| Part # | Description | Shipping Weight |
|------------|---------------------------------|-----------------|
| 7901243000 | Foot Ring for Leg Joint Comfort | 6 lbs. |





Exam Stools

Construction Specifications

Construction Specifications

| | |
|---------------------|--|
| Standard Dimensions | Seats: 15" diameter Polyfoam padding: 3" thick Height Range: 18" – 33" |
| Materials | Bases: Aluminum or Black Composite Seats: 74 colors of upholstered Vinyl covering polyfoam padding. Clifton/Passaic Stool: 304 stainless steel |
| Gauge | Clifton Stool Legs: 18GA; Clifton Stool Seat: 14 GA Passaic Stool Legs : 18 GA |
| Welding | Clifton/Passaic Stool: All welding by heliarc method and exposed welds are ground smooth. |
| Finish | Clifton/Passaic Stool: Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish. |
| Casters | 2" hard rubber casters or 2" twin wheel plastic casters. |
| Notes | MR conditional stool available. See MR conditional catalog. |

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

THE BLICKMAN DIFFERENCE


Vinyl Colors

- We offer an array of vinyl colors to fit your needs. Pick from 74 different colors that range from neutrals to cool tones to warm tones.
- Use 100% vinyl face for exam/patient rooms to ensure comfortable and easy-to-clean seating surfaces.
- Available for use in sterile environments. Antibacterial, Antimicrobial, and Mildew Resistant. Bleach cleaning allowed.

Key Product Features



Vinyl Colors

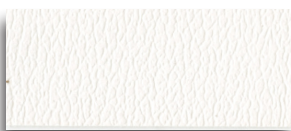

Snow White - 21



Platinum - 83




Parchment - 44




Coffee Cream - 81



Chocolate - 05

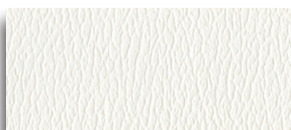

Super White - 06

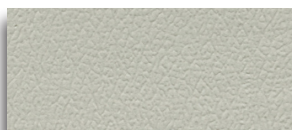


Alabaster - 38



Camel - 66



Sand - 09

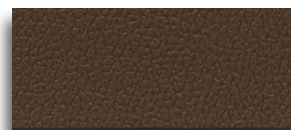


Chestnut - 61



Adobe - 48




Dove - 60

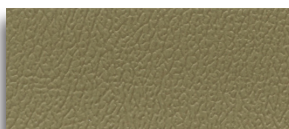


French Brown - 13




Light Maple - 14

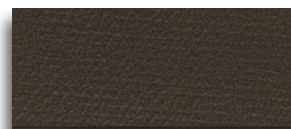

Mocha - 73

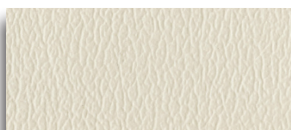

Moonstone - 16

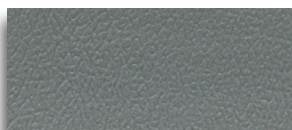


Franciscan Gray - 85



Taupe - 18



Brown Sugar - 19

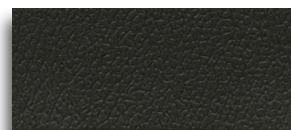

Espresso - 20



Quail - 87



Gunmetal - 02



Eucalyptus - 84




Cinnamon - 23



Black - 25


Vinyl Colors

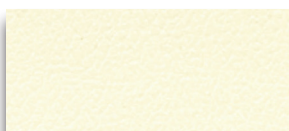

Summer Haze - 26



Light Lime - 27

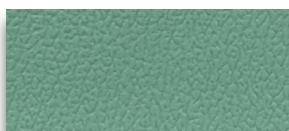

Celadon Green - 28



China Green - 29




Cosmos Blue - 30




Corn Silk - 31



Greenery - 32

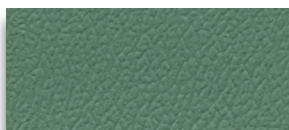


Jade - 33

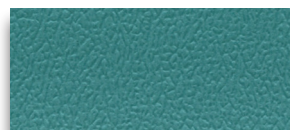

Turquoise - 34



Space Blue - 35

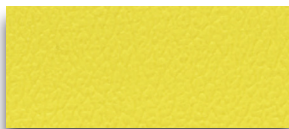


Lemon Peel - 36




Leaf Green - 37

BLICKMAN

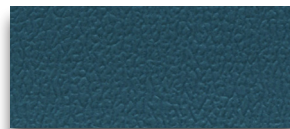


Spruce - 39

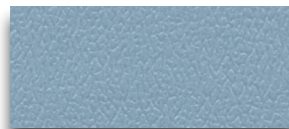


Aero Blue - 40




Sunflower - 41

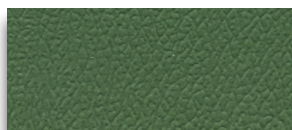


Loden - 42

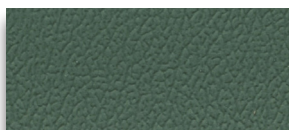


Viscaya Palm - 43

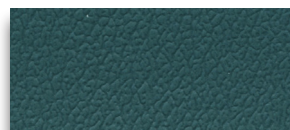


Grotto - 76

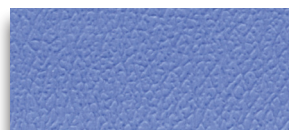


Azure Blue - 45

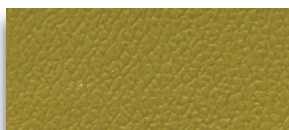

Tawny - 46

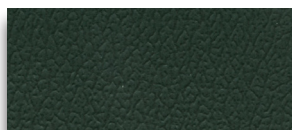

Yew Green - 47

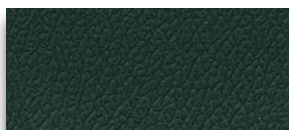

Forest Green - 77




Ocean Gray - 49




Regal Blue - 17


Vinyl Colors

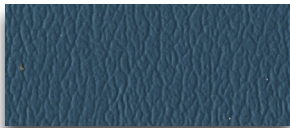

Colonial Blue - 51

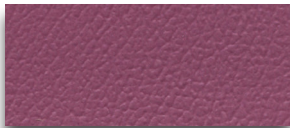


Crescent Cherry - 52

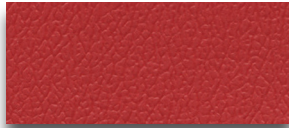

Tomato - 53




Dutch Orange - 54

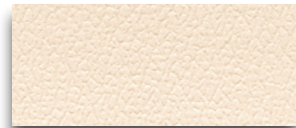


Shell - 55



Midnight Blue - 86

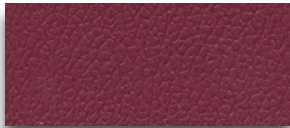

Raspberry - 57

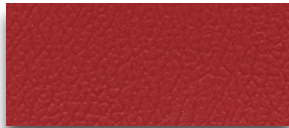

Lipstick - 58




Mandarin - 59




Clay - 82




Royal Blue - 80

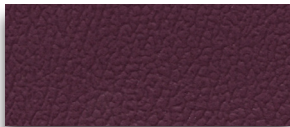

Grape - 56

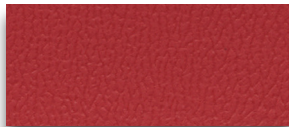

American Beauty - 63

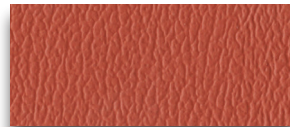

Copper - 64




Rosetan - 65




Academy Blue - 78



Dark Plum - 67

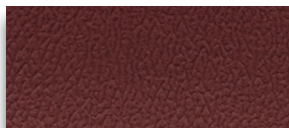


Bordeaux - 68

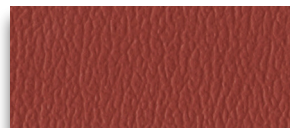

Cinnabar - 69




Petal - 70

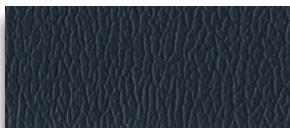


Navy Blue - 71




Deep Phoenician - 72

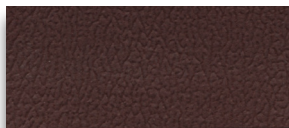

Burgundy - 15

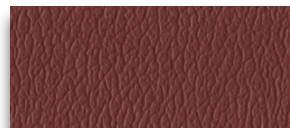

Paprika - 74

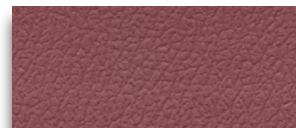

Tea Rose - 75



Vinyl Colors

Exam Stools

Cost:

- Stools come in all 74 colors at ***no extra cost.***

Composition:

- Antibacterial, Antimicrobial, and Mildew Resistant.
- **Face:** 100% Vinyl.
- **Backing:** 65% Polyester, 35% Rayon.

Specifications:

- **Weight:** 32.0 oz./lin. yd.

Abrasion Resistance:

- 100,000 Double Rubs.

Flammability:

- California Technical Bulletin 117-2013.
- FMVSS 302.
- NFPA 260 Class 1.
- UFAC Class 1.

Color:

- We provide an array of 74 Vinyl Colors.
- Dye lots may vary from samples shown.

Cleaning:

- Code: W
- Denim Dye Cleanable.
- Bleach Cleanable (80% Water, 20% Bleach).

Proposition 65 Warning: This product may contain chemicals known to the State of California to cause cancer, birth defects, and other reproductive harm.

*Final determination of the suitability of this product for application rests with the user.

BLICKMAN



OUR PROMISE

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality, and we pride ourselves on providing top notch and personalized service to ensure our customers are 100% satisfied. If you are ever not completely satisfied for any reason, we will make it right. And with our equipment backed by our Limited Lifetime Warranty, you can be certain you have made the very best choice for your healthcare facility.