

BLICKMAN



Scrub Sinks

Hands free multi-station operation

Includes:

- Windsor Scrub Sinks
- Lodi Scrub Sinks
- Accessories

OUR MISSION

Customer Focus

Exceed customer expectations through a commitment to continuous improvement of our products and services.

Quality Products

Offer the marketplace a broad range of high quality, cost effective, and durable products of superior craftsmanship.

Organizational Excellence

Strive for excellence in all we do by fostering a participative, team based environment, which nurtures the active involvement of all members of the organization.

Marketplace Recognition

Earn the respect, admiration, and loyalty of our people, customers, and suppliers by treating them with dignity and respect, embracing diversity, and developing and maintaining long term partnerships.

Collaborative Workplace

Encourage all employees to work safely and efficiently, contribute innovative ideas, and personally benefit from the ongoing success of the company.


The logo for Blickman, featuring the word "BLICKMAN" in a bold, white, sans-serif font. The text is contained within a dark blue rectangular box with a thin white border. The box has a slight 3D effect with a lighter blue shadow on the right side.

THE BLICKMAN DIFFERENCE

Scrub Sinks

- Our two scrub sink models, the Windsor and the Lodi, are each available with single, dual, and triple stations. A four station unit is also available in the Windsor model.
- They have deep sink bowls with sound deadened enclosed double walled panels that are angled in to divert water splashes along with spouts positioned toward the rear of the sink.
- All operation options are hands free. The Lodi scrub sinks have infrared water control and knee action water and soap dispensers, and the Windsor scrub sinks come standard with knee action control, but they are also available in foot, elbow, and infrared sensor controls.

Key Product Features

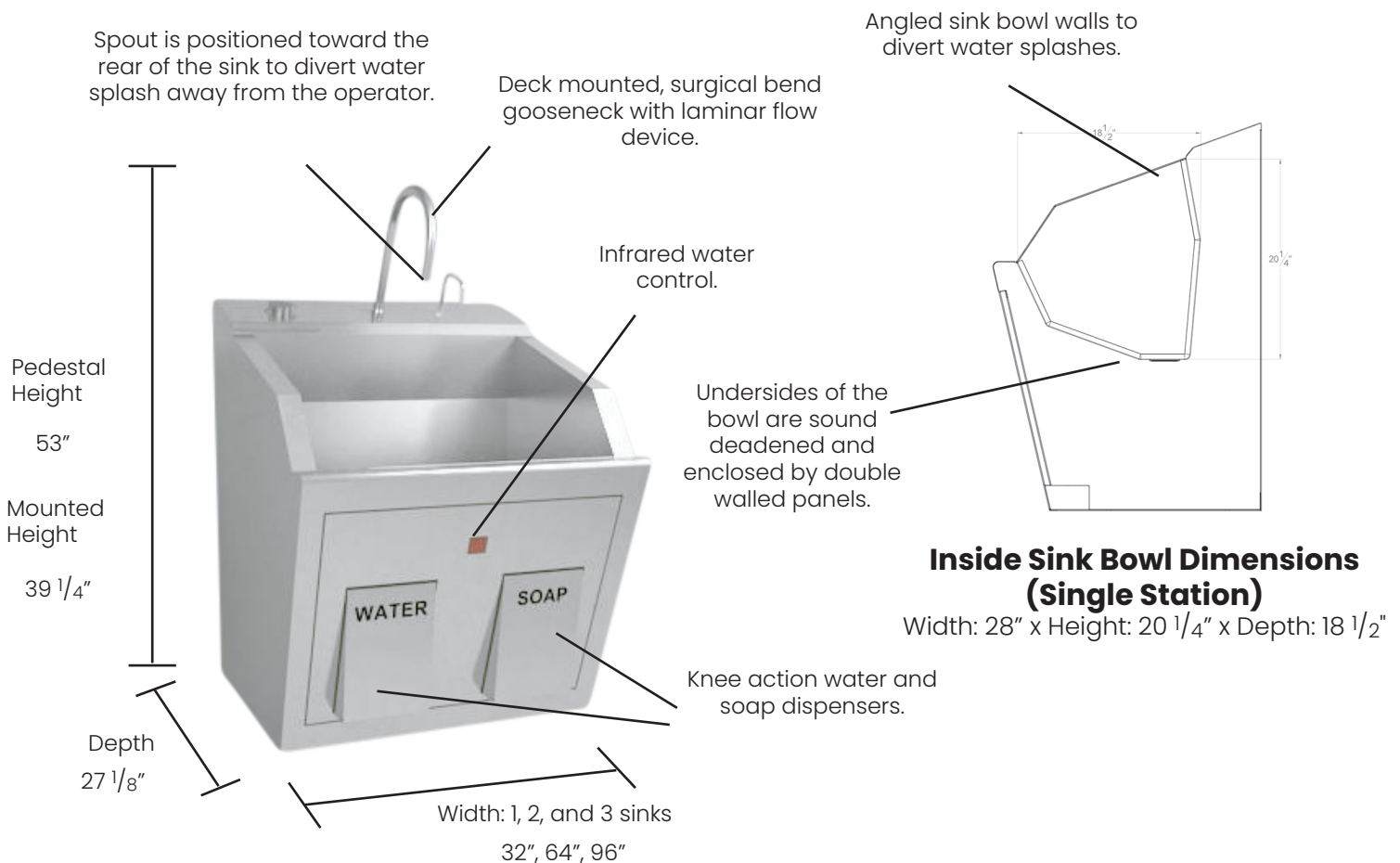

Scrub Sinks

Windsor

WINDSOR SCRUB SINKS

Description: Our Windsor scrub sinks are all wall mounted and utilize one sink bowl, undivided. Choose from our four different hands free operations: knee, foot, elbow, or infrared.

Inside Sink Bowl Dimensions

Width: 21 3/4" x Height: 8" x Depth: 14 1/2"



Width
25"



Width
47"



Width
69"



Width
91"

Name	1 Station	2 Stations	3 Stations	4 Stations
Knee Action				
Model #	7878SS-KA	7879SS-KA	7880SS-KA	7881SS-KA
Item #	1317878000	1317879000	1317880000	1317881000
Shipping Weight	65 lbs.	145 lbs.	190 lbs.	305 lbs.
Foot Action				
Model #	7878SS-FO	7879SS-FO	7880SS-FO	7881SS-FO
Item #	1317878003	1317879003	1317880003	1317881003
Shipping Weight	65 lbs.	145 lbs.	190 lbs.	305 lbs.
Elbow Action				
Model #	7878SS-EA	7879SS-EA	7880SS-EA	7881SS-EA
Item #	1317878004	1317879004	1317880004	1317881004
Shipping Weight	65 lbs.	145 lbs.	190 lbs.	305 lbs.
Infrared				
Model #	7878SS-IF	7879SS-IF	7880SS-IF	7881SS-IF
Item #	1317878001	1317879001	1317880001	1317881001
Shipping Weight	65 lbs.	145 lbs.	190 lbs.	305 lbs.



Scrub Sinks

Lodi

LODI SCRUB SINKS

Description: Our Lodi scrub sinks come standard with both infrared water control and knee action water and soap dispensers. All units are available in either pedestal or wall mounted configurations.



Width
32"



Width
64"



Width
96"



Height
53"



Height
39 1/4"

Name	1 Station	2 Stations	3 Stations
Pedestal			
Model #	9881SSP	9882SSP	9883SSP
Item #	1339881P00	1339882P00	1339883P00
Shipping Weight	275 lbs.	420 lbs.	600 lbs.
Wall Mounted			
Model #	9881SSW	9882SSW	9883SSW
Item #	1339881W00	1339882W00	1339883W00
Shipping Weight	215 lbs.	360 lbs.	570 lbs.



Scrub Sinks

Accessories

LODI SCRUB SINK ACCESSORIES

- Factory Installed Options -



Digital Timer



Eyewash

- For Wall Mount Installation -



Chair Carrier
For Wall Mount Installation

# of Stations	One	Two	Three
Part #	4301224751	4301224752	4301224753
Shipping Weight	47 lbs.	119 lbs.	141 lbs.

Scrub Sinks

Construction Specifications

Construction Specifications

Standard Dimensions	Lodi: (W) 32", 64", 96" x (D) 27 1/4" x (H): 39 1/4", 53" Internal Sink (Single Station): 28" x 18 1/2" x 20 1/4" Windsor: (W) 25", 47", 69", 91" x (D) 18 1/2" x (H): 15 3/4" Internal Sink: 21 3/4" x 8" x 14 1/2"
Materials/ Gauge	304 stainless steel Lodi: Mix of 16GA and 14GA Windsor: 16 GA Double/Triple Station Lodi Scrub Sinks: Plastic Dividers
Welding	All welding by heliarc method and exposed welds are ground smooth.
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish.
Electrical	Available in standard 100/120 VAC input or 220/240 VAC. Units are FDA registered and conform to standards for Safety Medical Electrical Equipment UL 606011 (US) and CSA C22.2 #601.1M90 (Canada). Units are cord connected with a hospital grade cap.
Casters	N/A

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.



OUR PROMISE

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality, and we pride ourselves on providing top notch and personalized service to ensure our customers are 100% satisfied. If you are ever not completely satisfied for any reason, we will make it right. And with our equipment backed by our Limited Lifetime Warranty, you can be certain you have made the very best choice for your healthcare facility.