Designer Series Fluid Warmers

The patient-friendly Designer Series allows you to put aesthetically pleasing warming units right where they are needed, allowing nurses and aides to maximize their time. The lighting inside the warmers softens the patient's room, placing more emphasis on comfort and care.

Safety is a priority with Enthermics Designer Series fluid warmers. Ensure fluids are safely warmed to recommended temperatures through the WarmSafe™ zone heating system. WarmRight® technology provides the right temperature for the right product to meet healthcare guidelines for IV and irrigation fluids.

WarmSafe Zone Heating: Heat is evenly balanced throughout the chamber to provide safe fluid warming.

Convenient: The hands-free, push button door is designed to make loading easier when carrying an armful of fluids.

Service: Our worldwide support team is dedicated to providing quality equipment. We walk with you every step of the way to maximize your return on investment.

Glass Door: Energy-efficient coated glass door reflects heat back into the cavity.

Made in the USA with a commitment to quality

An ISO 13485:2016 certified company

Standard Features

- · Windowed door
- Organizational basket
- Epoxy-coated exterior & interior surfaces
- Programmable control
- Interior lighting with a comforting blue glow
- Temperature management access point
- Four (4) non-skid rubber feet



28

DC400L:

28 (1) liter bags or 20 (1) liter bottles



Adjustable temperature range of

32°-40°C (90°-104°F) adjustable injection mode 32°-66°C (90°-150°F) adjustable irrigation mode

Accurate temperature

IRR +0/-1.67°C (-3°F)

INJ +0/-1.10°C (-2°F)





Configurations (select one)

Models

□ DC400L

Door Assembly

☐ Blue epoxy coated window door (standard)

Door Swing

☐ Right hinged (only)

Door Lock

☐ No lock (standard) ☐ Manual lock

Baskets

☐ 2 baskets (standard) ☐ 1 basket (optional)

Electrical

□ 120V 1PH □ 220V 1PH □ 230V 1PH

Accessories (select all that apply)

Bumper and Caster Assembly 3" (76mm) casters, 2 rigid, 2 swivel with Brake; overall height increases by 5" (127mm)

□ DC400L (5017547)

Caster Package with Plate 3" (76mm) casters, 2 rigid, 2 swivel with Brake; overall height increases by 5" (127mm)

□ DC400L (5017543)

Cart, Mobile

□ 3.5" (89mm) casters, 2 swivel; 2 swivel with brake (5024661); 36.5" (927mm) H x 26" (660mm) W x 26" (660mm) D

Anchorage (4 required)

☐ Anti-tip bracket (1014403)


Stacking Kits


It is recommended that fluid warming units be placed at or near the bottom of the stacked configuration whenever possible for maximum stability

□ DC400L stacked with DC400L or DC400 (5013024)





Specification

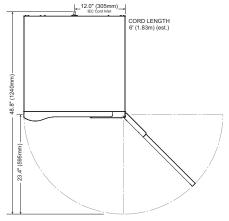

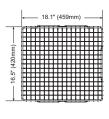



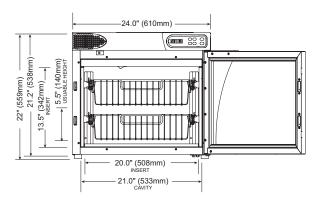

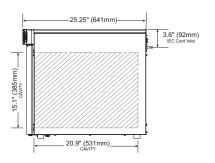



Top: 1"

Left: 1"

Right: 1" Back: 3"

DC400L

CLEARANCE



Exterior (H x W x D) Model

22"(559mm) x 24"(610mm) x 25.25 (641mm)

Interior (H x W x D) 15.25"(387mm) x 21"(533mm) x 21.5"(546mm) **Net Weight** 100 lbs (45 kg)

Ship Dimensions (L x W x H)*

31" x 31" x 32" 158 lbs DC400L

*Domestic ground shipping information. Contact factory for export weight and dimensions.



| | V | Ph | Hz | Α | kW | Plug Configuration |
|--------|-----|----|----|-----|-----|---|
| DC400L | 120 | 1 | 60 | 5.8 | 0.7 | NEMA 5-15P; 15A 125V Plug; Hospital Grade |
| | 220 | 1 | 50 | 2.7 | 0.6 | CEE 7/7 220-230V Plug |
| | 230 | 1 | 50 | 3 | 0.7 | BS 1363 Plug; CEE 7/7 220-230V Plug |

Ship Weight*



^{**}Electrical connections must meet all applicable federal, state and local codes.





IP-X0