

Be Independent.™



MOBILITY BEDRAIL™

BECAUSE A SWING OUT ARM MAKES STANDING EASIER THAN EVER



PRODUCT #5850



WATCH A DEMO





Includes a Bonus Organizer Pouch



Arm Pivots 180°



Adjustable Legs Provide Stability

FEATURES & BENEFITS

- **Swing Out Mobility Arm** – Pivots outward to provide improved walking stability.
- **Extendable Legs** – Provides added stability when standing up.
- **Ergonomic Cushion Handle** – Allows for easy transfer in and out of bed.
- **4-Pocket Organizer** – Provides storage space for handy items.
- **Anti-Slip Grips** – Secures rail in-between mattress and bed frame.
- **Height Adjustment** – Accommodates any home or hospital bed.
- **Easy Installation** – Installs in seconds with no tools required.

SPECIFICATIONS

Height Adjustment From Floor to Top of Handle: 21" - 44"
 Height Adjustment From Mattress Base to Top of Handle: 17.5" - 22.5"
 Rail Width: 19"
 Height Range From the Floor to Horizontal Support Structure Under Mattress: 13.5" - 21.5"
 Depth of Horizontal Support Structure Under Mattress: 25"
 Mobility Arm: swings out 19" away from bed
 Weight Capacity: 300 lbs
 Weight of Product: 10.5 lbs
 Made of: Steel frame with durable powder coat finish
 Package Dimensions: 16.25" x 24" x 2.875"
 Case of 4 Dimensions: 24.25" x 16.5" x 11.5"; 49 lbs


EASY ASSEMBLY (No Additional Tools Required)




Adjust Base Tubes to Accommodate Left or Right Side of Bed



Attach Base Tubes to Main Handle

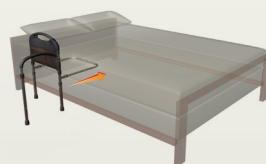

Insert Base Tubes Between Mattress and Box Spring

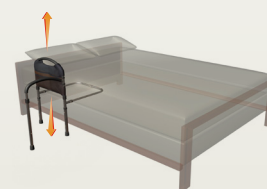

Adjust Height of Handle and Legs as Desired