





ODYSSEY

Extended Care Pressure Reduction Mattresses by Medline

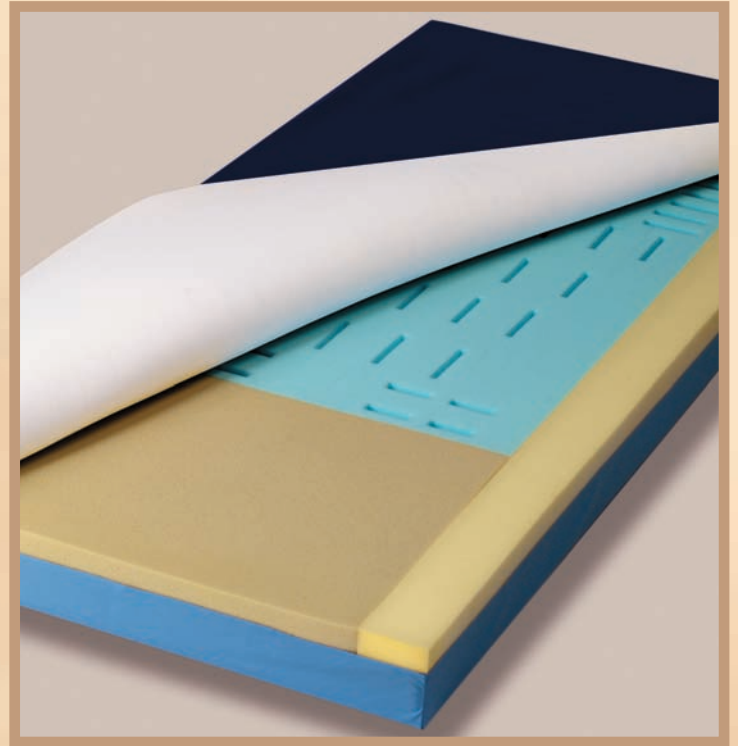

Engineered to minimize skin shear and pressure, while providing comfort and support for the entire body.

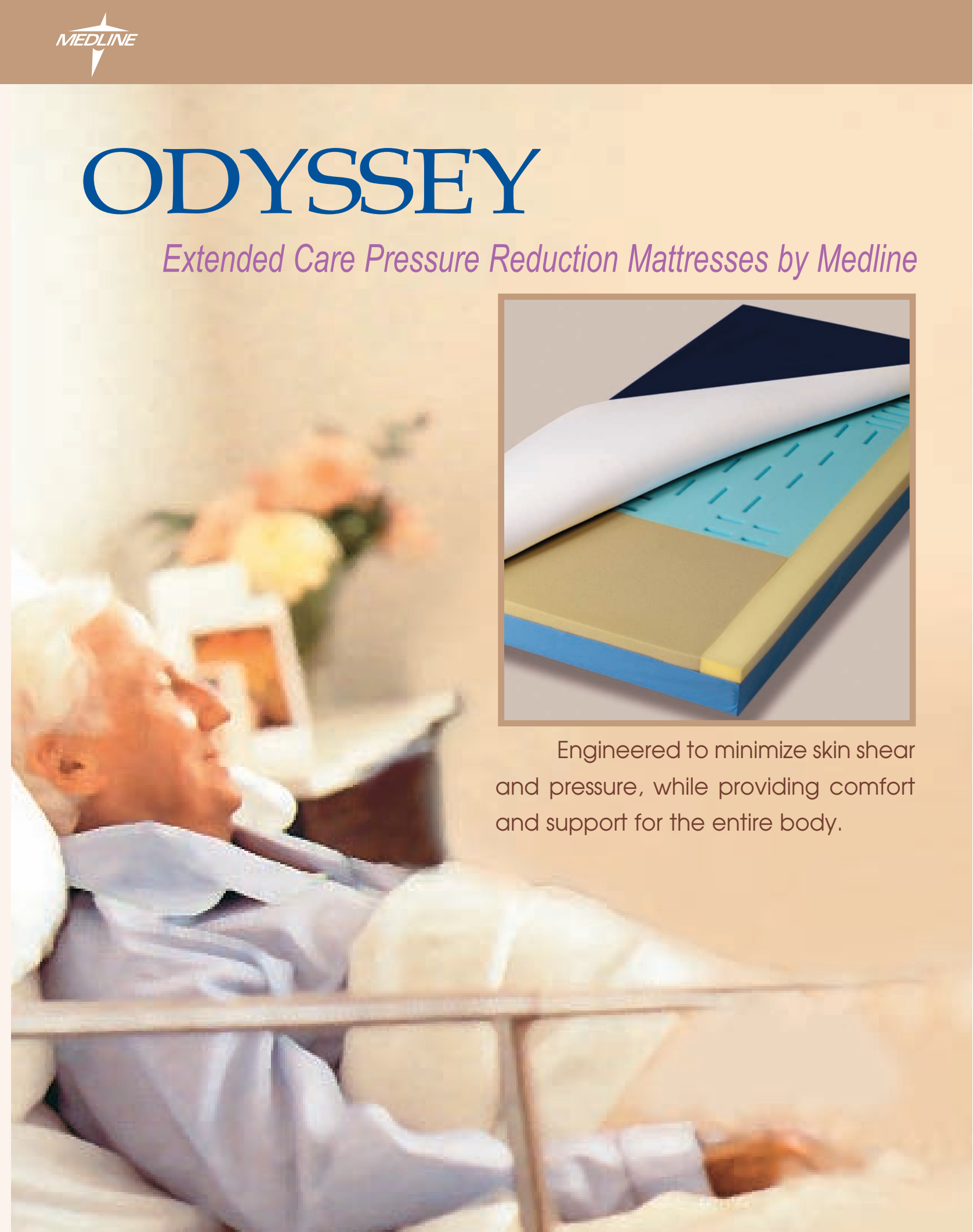

Odyssey Mattress Selection Guide

To assist you in choosing appropriate mattresses for your facility, Odyssey products are ranked below from greatest to lowest in therapeutic value.

Elite 5★

- Fluid-firm Visco foam top layer
- Double layer visco foam heel section
- Lowest-shear ProKNIT top cover
- Pressure Relief
- Helps eliminate powered mattress replacements. Use for stage 1-2 wounds, and stage 3 (without complications)

Zero-G 4★

- Nirvana Visco foam top layer
- Superior Pressure Reduction
- Helps eliminate the need for pressure reduction overlays. Use for stage 1-2 (without complications), and those with high risk of break-down

High-Res 3★

- High Resiliency foam top layer
- Moderate Pressure Reduction
- Helps eliminate overlays
- For moderate risk of skin breakdown

Extended Care 2★

- High Density foam top layer
- Basic Pressure Reduction
- Meets AHCPR guidelines for pressure reduction
- For low to moderate risk of skin breakdown

Nylox II Innerspring 2★

- Far more durable than vinyl mattresses, more comfortable and a ten-year limited Warranty
- For low to moderate risk of skin breakdown

Odyssey Innerspring 1★

- Economical general replacement



Odyssey Elite Mattress

5★ Our Highest Rated Mattress

Our NEW Odyssey Elite mattress features a Visco-Elastic conforming top layer which allows the patient to "sink" into the mattress and works to reduce pressure. Five anatomically-designed zones cradle bony prominences to lower skin shear and pressure. The ProKNIT top cover is fluid-proof, antimicrobial, and features all-way-stretch which helps prevent friction and skin shear. This mattress offers excellent durability for additional cost savings. Comes with a full 7 year warranty. Meets federal flammability requirements.



Item Number	Description
MDT360075 (A,B)	Odyssey Elite Mattress 36" x 75" x 6.5"
MDT360080 (A,B)	Odyssey Elite Mattress 36" x 80" x 6.5"
MDT360084 (A,B)	Odyssey Elite Mattress 36" x 84" x 6.5"

Pressure Reduction Mattresses from Medline

Odyssey Zero-G Mattress 4★ Performance

Our top of the line, the Zero-G, uses heat activated Nirvana foam coupled with a Nylex Ultra-Soft cover to conform to body contours. This helps reduce pressure on all areas of the resident's skin. The extremely high density of Nirvana foam makes it longer lasting. The Odyssey Zero-G is a great addition to your wound prevention program. Priced affordably, the Zero-G can help you reduce costs by eliminating the use of specialty overlays. Full five year warranty.

Item Number	Description
MDT230175 (A,B)	Odyssey Zero-G Mattress 36" x 75" x 6"
MDT230180 (A,B)	Odyssey Zero-G Mattress 36" x 80" x 6"
MDT230184 (A,B)	Odyssey Zero-G Mattress 36" x 84" x 6"

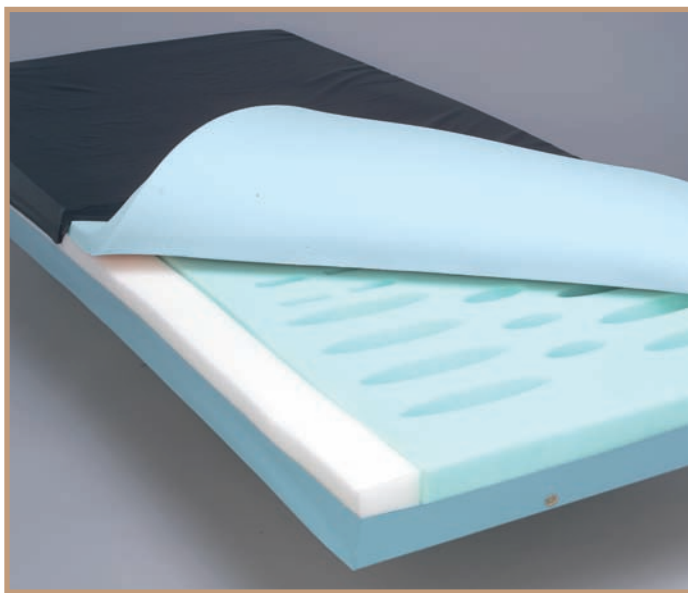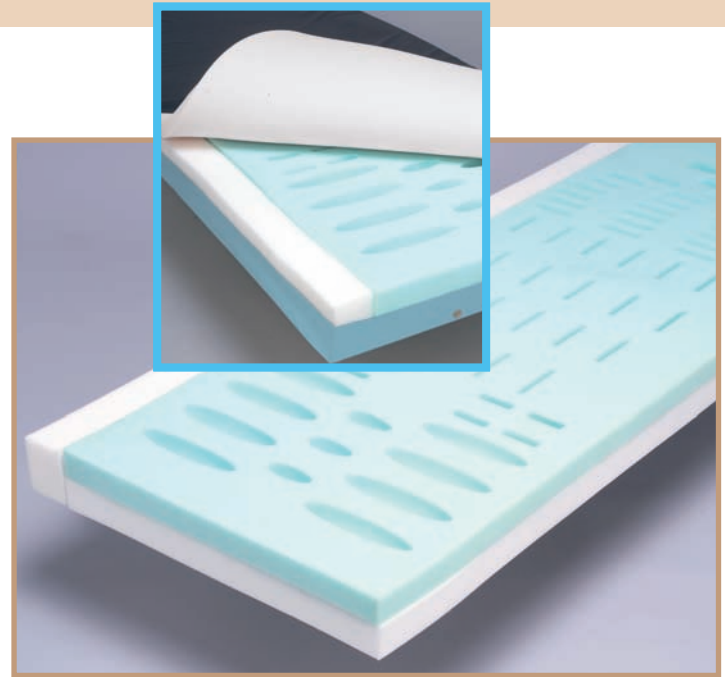

Odyssey High-Res Mattress 3★ Quality, Comfort and Pressure Reductions

This layered foam mattress combines a firm base support layer, extra firm edges and our exclusive Odyssey five-zone center foam layer. The top layer of foam is made of 1.5-inch thick high resiliency polyurethane foam and covered by our Nylex Ultra-Soft ticking. The High-Res is an excellent support surface for residents at risk for pressure ulcers. Full five year warranty.

Item Number	Description
MDT230275 (A,B)	Odyssey High-Res Mattress 36" x 75" x 6"
MDT230280 (A,B)	Odyssey High-Res Mattress 36" x 80" x 6"
MDT230284 (A,B)	Odyssey High-Res Mattress 36" x 84" x 6"

Medline Bumper Wedges & Gap Protectors

Item Number	Description		
MSC04200	Bumper Wedge for side rails	14" x 70"	2 ea/cs
MSC04205	Bumper Wedge for side rails	14" x 60"	2 ea/cs
MDT23NOGAP	Mattress/side rail gap protector	80" x 2"	1 ea/cs
MDT23NOGAP2	Mattress/head or footboard gap protector	36" x 2"	1 ea/cs
MDT23NOGAP3	Mattress/side rail gap protector	80" x 4"	1 ea/cs
MDT23NOGAP4	Mattress/head or footboard gap protector	36" x 4"	1 ea/cs

* See back cover for SAF-T-SIDE Mattresses





Odyssey Extended Care Mattress

2★ Comfort Plus

This layered foam mattress combines a firm base support layer, extra firm edges and our exclusive Odyssey five-zone center foam layer. The top layer of foam is 1-inch thick high-density polyurethane foam and covered by our Nylex Ultra-Soft ticking. Our most cost effective Odyssey mattress, this is your choice for super value. Full three year, plus two year prorated cover warranty.

Item Number	Description
MDT230375 (A,B)	Odyssey Extended Care 36" x 75" x 6"
MDT230380 (A,B)	Odyssey Extended Care 36" x 80" x 6"
MDT230384 (A,B)	Odyssey Extended Care 36" x 84" x 6"

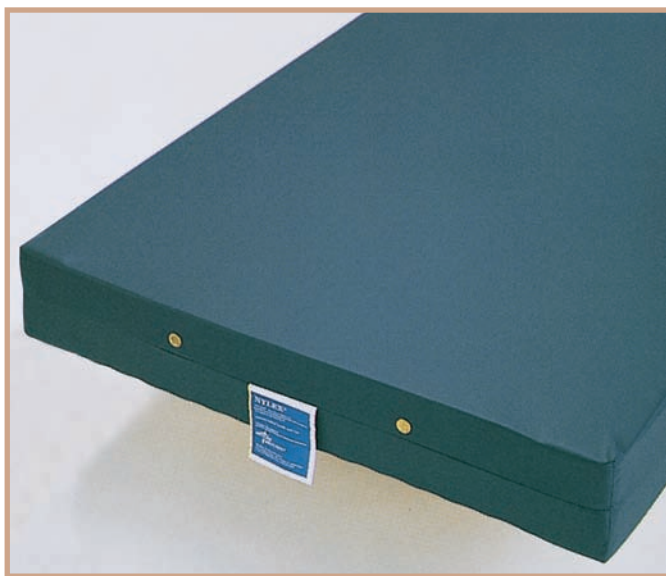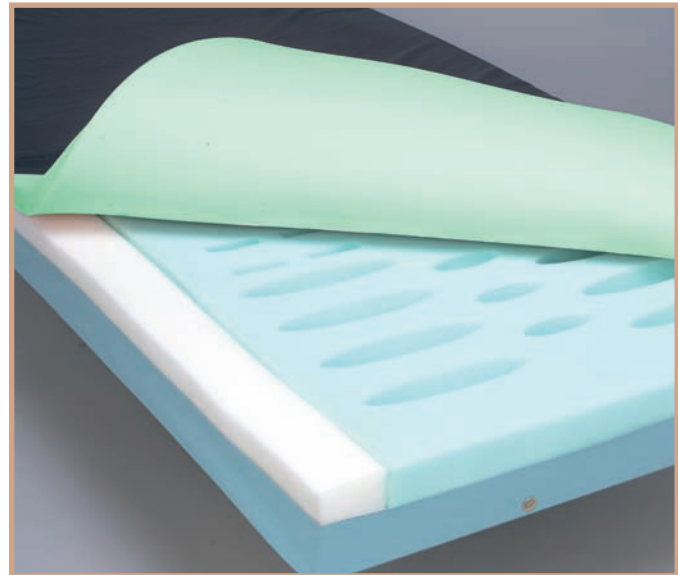

Odyssey Innerspring Mattress

1★ Best Value

Made with an inner core of high-density foam on both sides of a tempered steel, low profile innerspring unit. This mattress provides the support you want with the comfort you need in long term care. Made with firm edges, this mattress is turnable and flipable to promote long mattress life. The comfort designed Nylex Ultra-Soft cover is made with an inverted side seam for easy cleaning.


Item Number	Description
MDT230475 (A,B)	Odyssey Innerspring Mattress 36" x 75" x 5.5"
MDT230480 (A,B)	Odyssey Innerspring Mattress 36" x 80" x 5.5"
MDT230484 (A,B)	Odyssey Innerspring Mattress 36" x 84" x 5.5"

Nylex II Innerspring Mattress

2★ Our Top of the Line Innerspring

Built to last, with convoluted foam around a 308-coil, low-profile steel innerspring unit, the Nylex II provides excellent support and comfort. Medline's Nylex ticking is fluid-proof, antimicrobial, flame-retardant and stain-resistant. It's easy to understand why the Nylex II is our most popular pressure reduction replacement mattress.

Item Number	Description
MDT231275 (A,B)	Nylex II Innerspring Mattress 36" x 75" x 6"
MDT231280 (A,B)	Nylex II Innerspring Mattress 36" x 80" x 6"
MDT231284 (A,B)	Nylex II Innerspring Mattress 36" x 84" x 6"

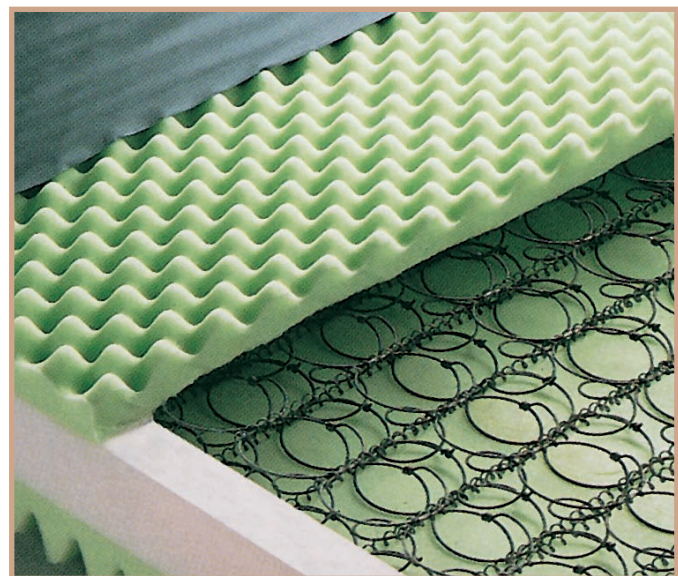

When ordering any of the mattresses featured in this brochure, please add the appropriate letter code after the item number to indicate the total quantity ordered:

Quantity	Code
1-5	A
6+	B

For example: MDT230180A would indicate an order of 1 to 5 each 36" x 80" Odyssey Zero-G Mattresses.

** To add Saf-T-Sides (see back cover) use code REA for 1-5 and REB for 6 or more. i.e., MDT230180REB would be 6+ Odyssey Zero G's with Saf-T-Sides

As the leading supplier of Pressure Reduction mattresses to long-term care facilities, Medline is proud to introduce our latest improvement, the Odyssey series of mattresses. As was the case sixteen years ago with our introduction of fluid-proof Nylex ticking and multi-layered antimicrobial foam, Odyssey sets the new standard for health care mattresses.

Our customers trust us to help them achieve their wound prevention and treatment goals. We have earned this trust by providing excellent quality, affordable mattresses, and the expertise to help determine the appropriate model for every care situation.

Due to infection control requirements of the health care market, we design all of our mattresses to be fully antimicrobial. All of our ticking and foam is treated, inside and out, with Ultra-Fresh™ antimicrobial agent. This treatment, effective for the life of the mattress, provides protection against a wide spectrum of organisms including MRSA, and the odors that may be associated with bacteria.

Nylex Ultra-Soft, our most supple cover material yet, goes even further to help reduce friction and shear. Nylex Ultra-Soft provides exceptional resident comfort from its softness and the elimination of heat buildup. It also offers excellent stain resistance and ease of cleaning for your staff.

You can be sure Medline will continue to give you the most health care mattress for your money!

DESIGNED FOR COMFORT AND PERFORMANCE

- Each of the Odyssey mattresses is designed with an advanced center layer. Also, five different zones of anatomically designed dye cuts help minimize pressure and skin shear from head to foot.



- Larger cuts at the foot end cradle the fragile bony prominences of heels and ankles. Smaller cuts in the trunk help lower shearing forces while providing skeletal support.

- The advanced design helps the Odyssey line promote good posture while relieving pressure and providing comfort.



- Inverted side seam seals top and bottom cover materials to provide superior protection against fluid penetration. Brass vents allow the mattress to breathe. Slip-resistant vinyl bottom material withstands day-in, day-out institutional use.



For More Information
Call Toll-Free:
1-800-MEDLINE
(1-800-633-5463)



Medline Extended Care Mattress Options & Accessories



■ Removable zippered cover with barrier flap to prevent fluid penetration is available on all layered foam mattresses. Ask to see our optional ProKNIT all-way-stretch top cover. This knitted synthetic fabric with a fluid-proof coating is soft, smooth and built to last. ProKNIT is antimicrobial, flame-retardant, and stain-resistant.



■ Four evacuation handles can be added to any of the mattresses in this brochure. Made of high-strength nylon webbing and box-tacked to the bottom of the mattress cover, these handles make it easier to transport the mattress.



■ Saf-T-Side raised edges are designed for use when bed side rails are not raised. Featuring an easy exit-entry cutout which enables residents to sit up and get in and out of bed when they want to, our firm Saf-T-Side raised edges help to prevent resident falls. Available on all layered foam mattresses.



■ All Medline mattresses are flame retardant*. But any of our mattresses can provide an even higher level of safety with the addition of our optional fire barrier. Fireproof netting surrounds the foam core of the mattress and meets stringent open-flame-testing standards including California Technical Bulletin #129 and City of Boston Fire Code IX-11.

*Flame Retardant Information

All foam components used in Medline mattresses featured in this brochure meet or exceed California Technical Bulletin #117. All mattress covers meet or exceed NFPA 701 (Small Scale Test). In addition, all mattresses featured pass the Cigarette Test (CCF-4-72). With the addition of the optional fire barrier, Medline mattresses meet California Technical Bulletin #129 and City of Boston Fire Code IX-11.

Warranty

Medline Mattresses are Warranted by the Manufacturer as follows:

Under Medline's full warranty, the mattresses will be repaired at the manufacturer's sole option. Additionally, no returns or replacements will be allowed without authorization from Medline Industries, Inc. Upon authorization, prepaid freight will be the responsibility of the end user. The foregoing Warranties of the manufacturer are in lieu of all other warranties either expressed or implied, including, without limitation, any implied warranties, merchantability or fitness for a particular purpose. The manufacturer is not liable for consequential or incidental damage due in whole or in part to the use of its products except as otherwise agreed or as required by law. Please refer to individual product listings for length of warranty.



Medline Industries, Inc.

One Medline Place

Mundelein, IL 60060.4486

Toll-Free: 1-800-MEDLINE

www.medline.com