


Protekt™ 300

Pressure Redistribution Foam Mattress

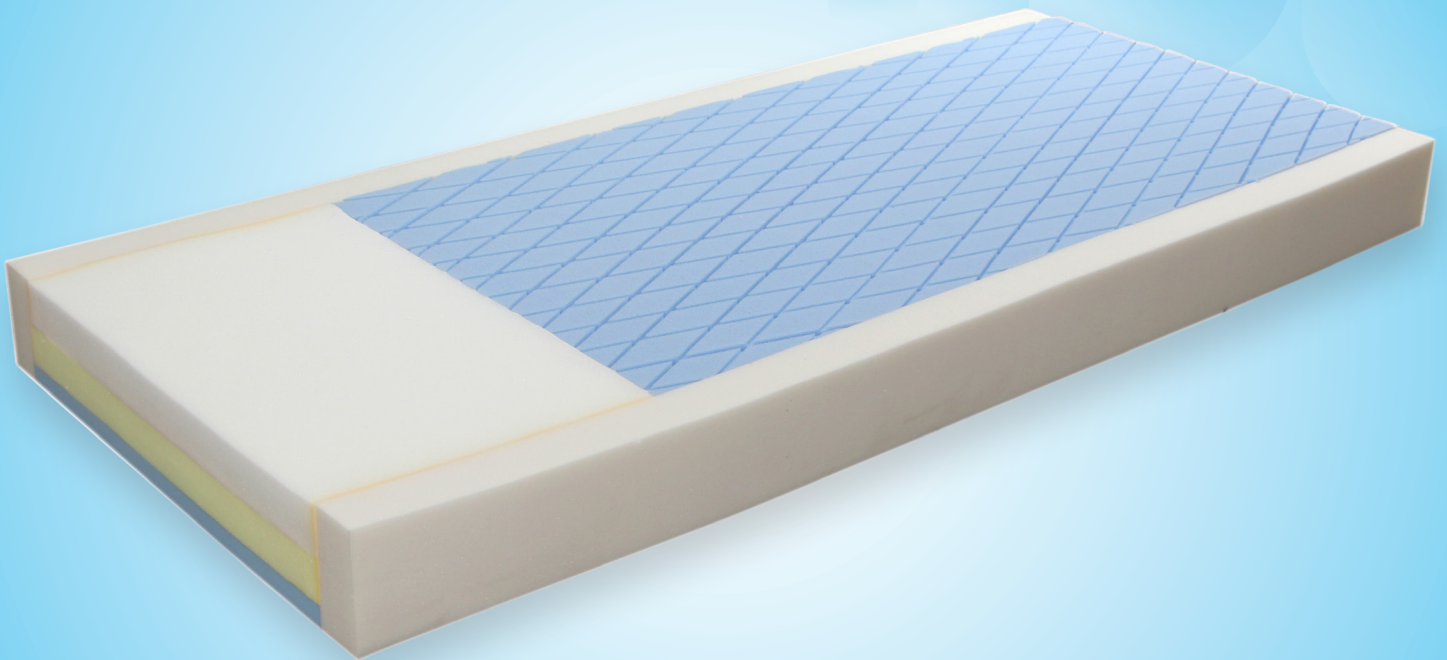

 **Proactive**
medical products

Medicare
Code
E0184

270 Washington Street Mount Vernon, NY 10553

Tel: 855-BE-PROACTIVE (855-237-7622) • Fax: 845-205-6005

www.proactivemedical.com • sales@proactivemedical.com

The Protekt™ 300 provides exceptional pressure redistribution. This mattress is ideal for moderate to high risk patients with close attention paid to the vulnerable sacral and heel area by utilizing ultra high density foam.

FEATURES AND BENEFITS:

- Multi-tiered pressure redistribution mattress with top layer of die-cut high density foam in the head and torso section and an exceptionally high resilient foam in heel section provides an unparalleled pressure redistribution support surface.
- Middle layer of heel section has air channel cut foam which provides additional pressure redistribution in heel area and air circulation.
- PTSS Rail System provides additional support, allowing easy patient transfer and edge of bed sitting.
- Top cover is a 4-way stretch nylon which reduces perspiration and provides a low shear friction surface that is anti-microbial and water resistant.
- Anti-Contamination Flap covers the zipper reducing the possibility of fluids contaminating the foam core.
- Designated head and foot sections.
- Fire sleeve to meet Federal Fire Code 16 CFR part 1633.
- Weight Capacity: 350 Lbs.
- Warranty: 7 year, Non Prorated Warranty.

81031	Protekt 300	Pressure Relieving Foam Mattress	36"x76"x6"
81032	Protekt 300	Pressure Relieving Foam Mattress	36"x80"x6"
81033	Protekt 300	Pressure Relieving Foam Mattress	36"x84"x6"
81034	Protekt 300 Plus	Pressure Relieving Foam Mattress w/3" Raised Rails	36"x76"x9"
81035	Protekt 300 Plus	Pressure Relieving Foam Mattress w/3" Raised Rails	36"x80"x9"
81036	Protekt 300 Plus	Pressure Relieving Foam Mattress w/3" Raised Rails	36"x84"x9"



Optional raised side rails



270 Washington Street Mount Vernon, NY 10553

Tel: 855-BE-PROACTIVE (855-237-7622) • Fax: 845-205-6005

www.proactivemedical.com • sales@proactivemedical.com