

Protekt™ Aire 6000

Low Air Loss/Alternating Pressure Mattress System



 **Proactive**
medical products

Medicare
Code
EO277

270 Washington Street Mount Vernon, NY 10553

Tel: 855-BE-PROACTIVE (855-237-7622) • Fax: 845-205-6005

www.proactivemedical.com • sales@proactivemedical.com

The Protekt™ Aire 6000 is the preferred choice by clinicians treating stage I-IV pressure ulcers. The combination of low air loss and alternating pressure therapy's provide an unparalleled pressure redistribution surface.

FEATURES AND BENEFITS

MATTRESS & COVER

- Unique cell-on-cell design ensures that the mattress will hold air in the lower cells in the event of a power failure.
- Transport connector locks in air for easy patient transport.
- Quilted top cover provides comfortable sleeping surface.
- Bottom cover has 3 non-skid rubber panels to prevent mattress from slipping.
- Dimensions: 35"x80"x8".

PROTEKT™ AIRE 6000 PUMP

- Powerful yet quiet hi-flow 8 LPM (liters per minute) pump.
- Dual therapies; low air loss and alternating pressure.
- 4 alternating times to customize therapy.
- Static mode provides firm support for egress or ingress from mattress or side edge seating.
- Seat-Inflate feature helps protect the patient when they are in the fowler position.
- Audible and visual alarms in case of low pressure.
- Auto firm mode for quick installation.
- UL/CSA approved power cord.
- Weight Capacity: 350 lbs.
- 18 month warranty.



Transport Handles




Easy Valve Exchange



Low Air Loss



Transport Connector

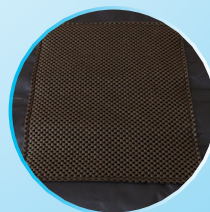

Non-skid Bottom



CPR Release



Quilted Cover

80060	Protekt™ Aire 6000	8" Alternating/Low Air Loss Mattress System
80061	Protekt™ Aire 6000 Pump	Pump only for Protekt™ Aire 6000
80062	Protekt™ Aire 6000 Mattress	Mattress only for Protekt™ Aire 6000
80063	Protekt™ Aire 6000 Cover	Cover only for Protekt™ Aire 6000
80064	Protekt™ Aire 6000 Bladder	Individual Bladder for Protekt™ Aire 6000
80065	Protekt™ Aire 6000 Raised Rails	8" Alternating/Low Air Loss Mattress System w/Raised Rails
80066	Protekt™ Aire 6000 84"	8" Alternating/Low Air Loss Mattress System-84"
80067	Protekt™ Aire 6000 84" Raised Rails	8" Alternating/Low Air Loss Mattress System-84" w/Raised Rails
80068	Protekt™ Aire 6000 84" Mattress	84" Mattress only for Protekt™ Aire 6000
80069	Protekt™ Aire 6000 84" Cover	84" Cover only for Protekt™ Aire 6000



270 Washington Street Mount Vernon, NY 10553

Tel: 855-BE-PROACTIVE (855-237-7622) • Fax: 845-205-6005

www.proactivemedical.com • sales@proactivemedical.com